

Celexo

Hydro Gel



Celexo
Hydro Gel

Celexo

Skin Serum



Spec	3ml x 2vial / 1box
Type	Vial
Indication	Anti-inflammation, Skin barrier repair
Key Ingredients	Exo-Cica®, Hybrid HA (High&Low Molecular HA)

Celexo

Hydro Gel



Spec	2ml x 2syringe / 1box
Type	Syringe
Indication	Boost Collagen, Tissue Regeneration
Key Ingredients	Exo-Cica® High Molecular Weight HA

Celexo *Hydro Gel*



VOLUMIZING & HYDRATING SOLUTION

Tissue Regeneration and Hydration

GUARANTEED SAFETY

Clinically Validated Skin Regeneration

SKIN REGENERATION

Exo-Cica®



HYDRATION LOCK

HMW-HA



*High Molecular Weight Hyaluronic Acid

GUARANTEED SAFETY

Clinically Validated Skin Regeneration

Class II Medical Device

Certified safety

DNA-Free, Cell-Free Formula

No animal or
human-derived genetic material

Safe Application

Pre-filled syringe ensures sterile
and easy use

No Cold Chain Required

Stable at room temp,
with stands 121°C autoclave

Celexo Hydro Gel




2 mL x 2 syringe / 1box
[Approval No.] 25-4356 | [Storage] 1~30°C

KEY INGREDIENT **Exo-Cica**®


Next-generation bio-regenerative material

- **Cica leaves were harvested at Smart Farm** from ABio's own variety 'Hoin'
- **Extracted by ABio's own patented technology as ExoTraction**®
- **Cica exosome showed anti-inflammatory, cell regeneration, whitening and skin barrier improvement efficacy**

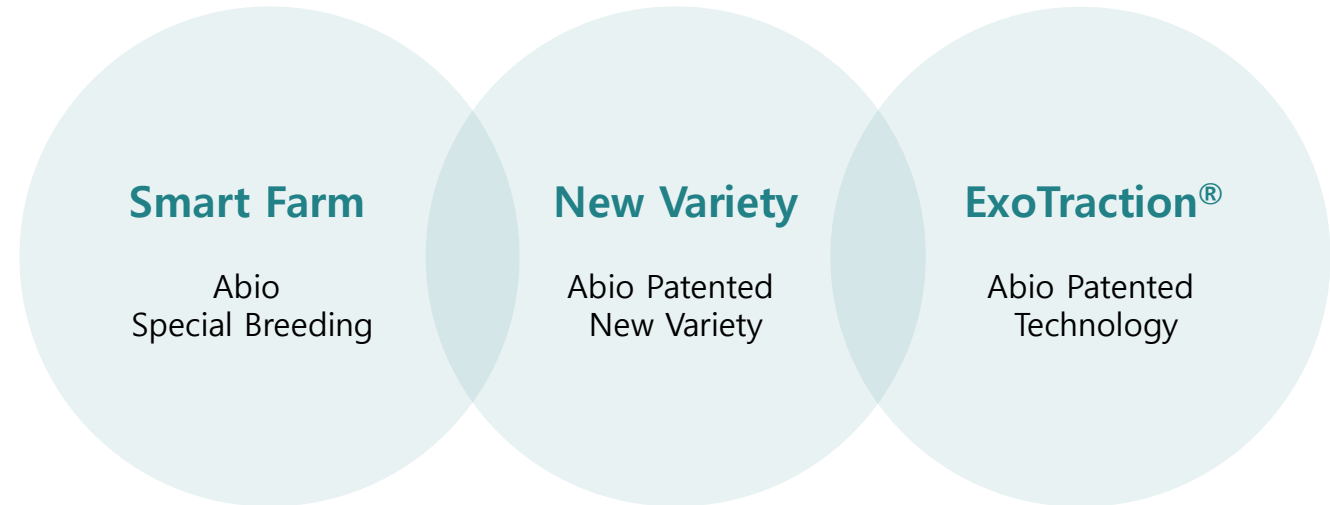


KEY INGREDIENT Exo-Cica®

Next-generation bio-regenerative material

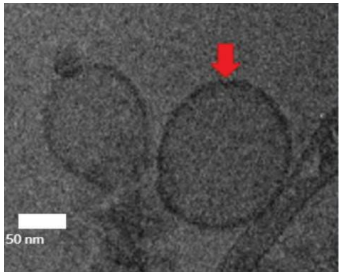

Cryo-TEM analysis image

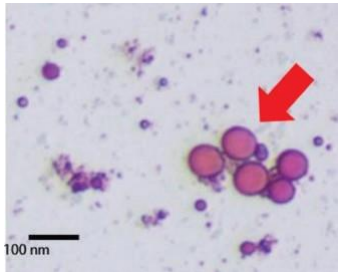

Dual vesicle structure through Oil Red O staining



A patented ingredient produced using A-bio Materials's proprietary **ExoTraction®** technology for exosome isolation and purification.

▪ miRNA Quantitative Analysis Result

Exo-Cica® miRNA analysis confirmed high-purity miRNA expression

Quantitative analysis of miRNA, **effective for regenerative and anti-inflammatory signaling**

Item	Result	Evaluation
miRNA concentration	110.3 ng/μl	Very high concentration(Excellent)
A260/A280	2.05	Excellent purity (minimal protein/nucleic acid concentration)
A260/A230	1.83	Excellent (low salt/organic compound contamination)

▪ Phosphatidic Acid (PA) Analysis

Item	Result
PA content	0.25%
Analysis Method	HPLC
Analysis Meaning	High PA content contributes to strengthening the skin barrier


Exo-Cica® **minimizes damage** by being cholesterol free

KEY INGREDIENT **Exo-Cica®**

Next-generation bio-regenerative material

■ Stability Evaluation

NTA analysis to confirm particle number and size distribution changes

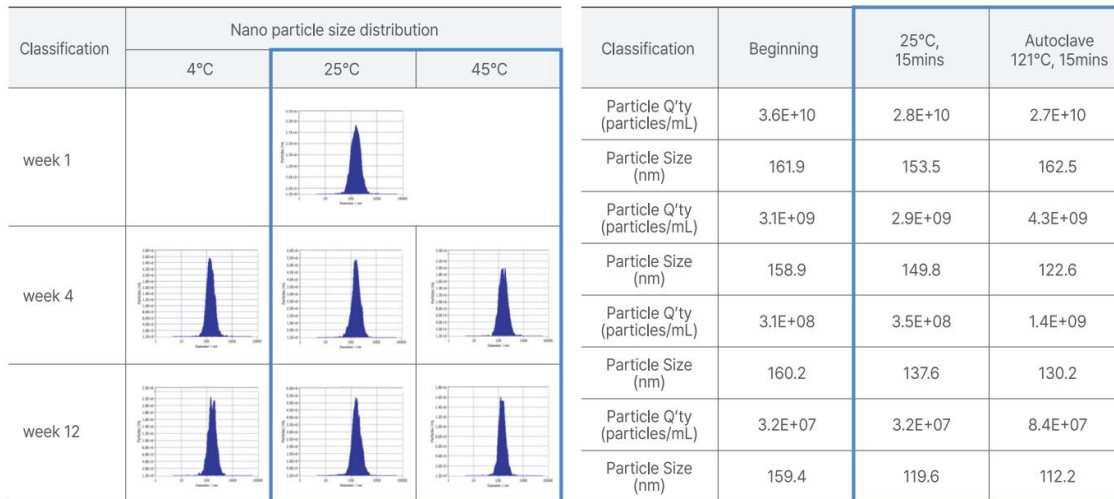


From week 0 to week 12, **remains stable up to 12 weeks** at various temperature while maintaining size and penetration.


■ Exosome Skin Penetration Test

Verification of penetration and cellular activation by Exo-Cica®

[NTA results]

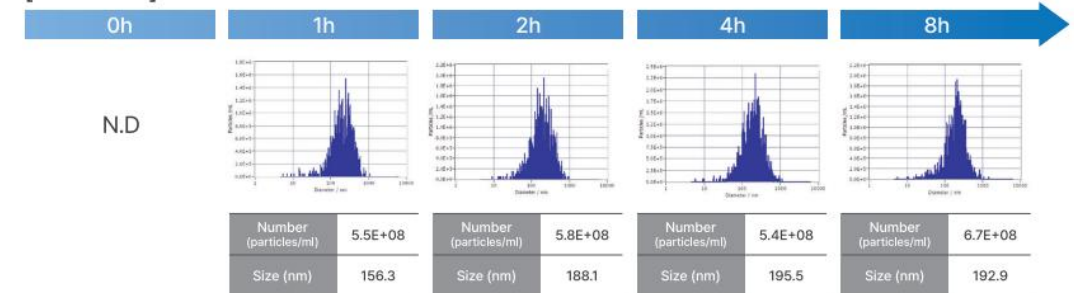


[Oil red O results]

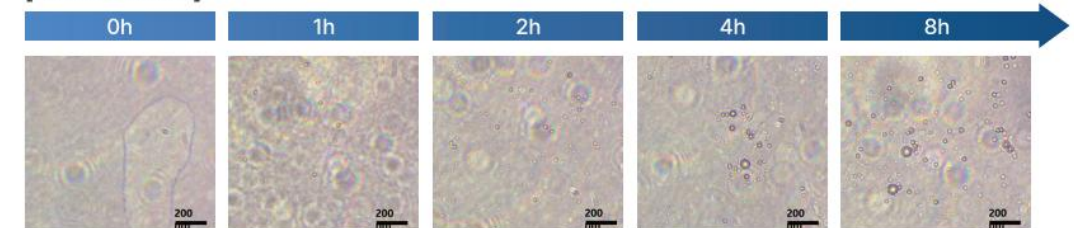

Microscopic images **confirm absorption** into skin cells over time
→ Proves that the active ingredients were delivered.

KEY INGREDIENT **Exo-Cica®**



Next-generation bio-regenerative material

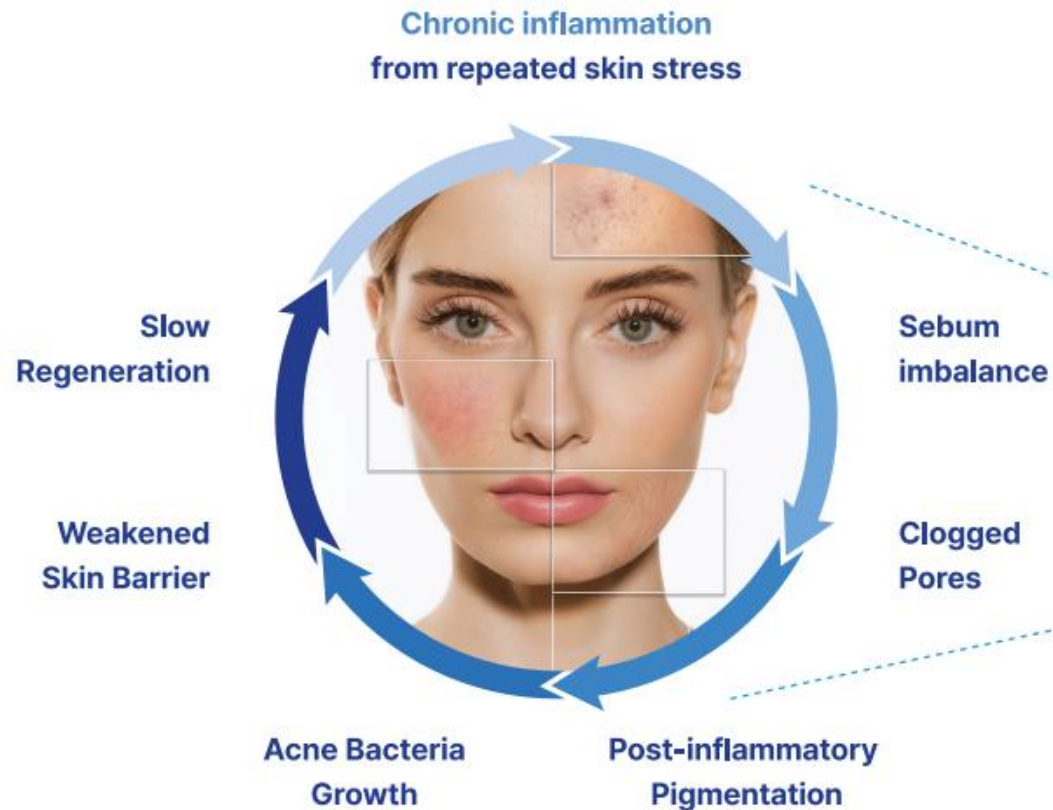

Skin Regeneration Solution via Anti-inflammatory Action

Cell Regeneration

Regenerative Solution **3-Stage Full Cycle of Regeneration**

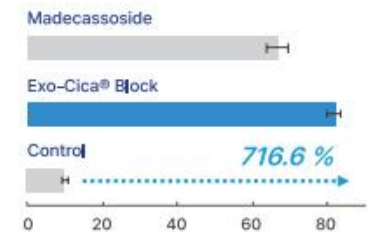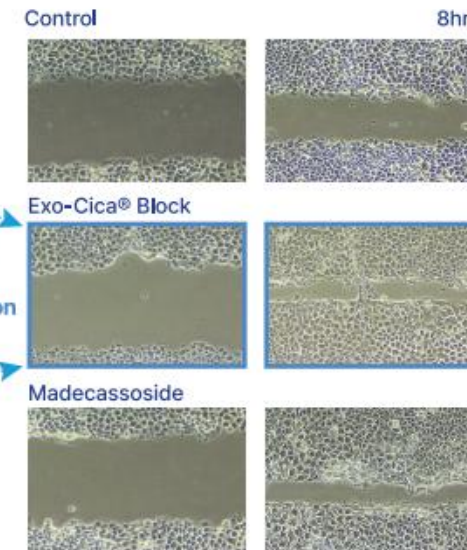


HDF Scratch Wound Recovery Test for Cell Regeneration Confirmation

KEY INGREDIENT Exo-Cica®


Next-generation bio-regenerative material

Anti-inflammation

Ability to **suppress inflammatory cytokine** site expression

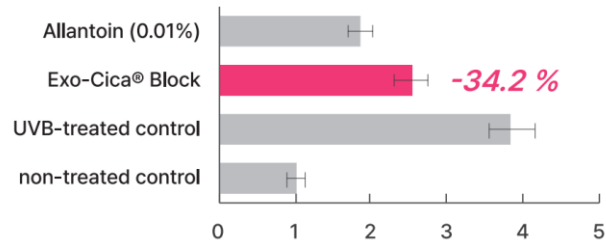


TNF-α mRNA expression

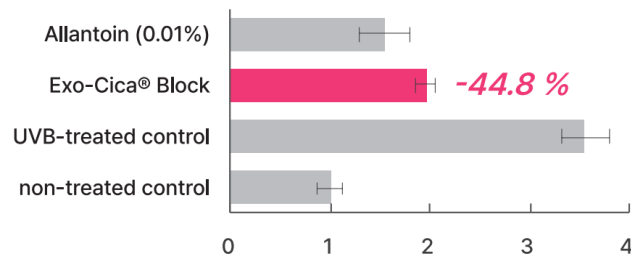


IL-6 mRNA expression

Collagen Synthesis & Melanin Suppression-Barrier Restoration


Complete Restoration: Collagen increase, Melanin suppression, Barrier restoration

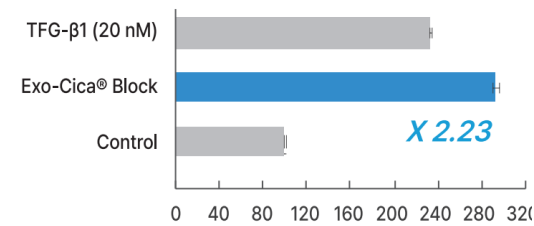


COL1A1 mRNA expression

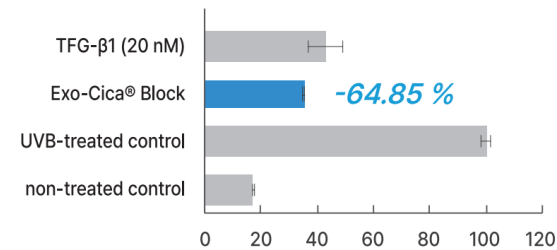

MMP-1 mRNA expression

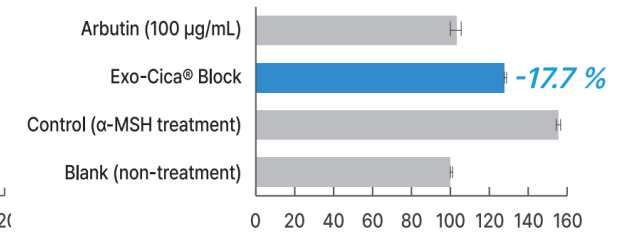

Melanin content (% of control)

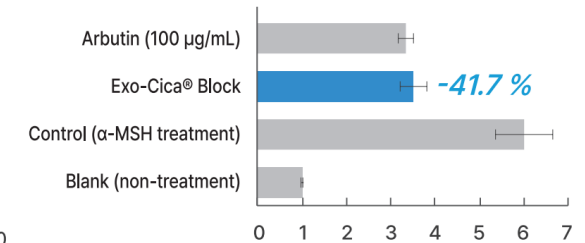

TRP-1 mRNA expression

KEY INGREDIENT Exo-Cica® vs. Cell culture media

Category	Centella Asiatica Extract	Callus Culture Extract	Exo-Cica®
Origin	Dried or fresh Centella(extract)	Plant cell callus	Fresch Centella leaf (New variety 'Hoin')
Extraction Method	Ethanol / Water extraction	Extraction after culturing	Exo-Traction® (cell membrane + ATPs)
Main Ingredient	Madecassoside, Asiaticoside, etx.	Antioxidants at basic levels	Exosomes + RAN + Proteins + Lipids
Exosome	X (Not an exosome)	X (Not an exosome)	High-purity exosome
Bioactive Compound Diversity	Medium	Low	Very High
Skin Penetration & Delivery	Medium	Low	Very High
Skin Regeneration	Medium	Medium	Very High
Key Differentiator	General skincare, widely used	Low-cost alternative, low efficacy	Patented exosome for skin booster – The next level evolution of Cica


General Cosmetic Ingredients vs. Skin Booster Grade High-Concentration Exosome

- **More than 2x response** to gene expression
- Activation of regenerative genes **by about 46%**
- **About 70% expression** of barrier genes
(e.g. FGF, LOX, H3S3T)
Reduction of aging-related genes **by average of 35%**


Celexo Hydro Gel forms a bioactive film that delivers exosomes and locks in moisture-accelerating recovery




In-Vitro Data | Skin Regeneration & Hydration

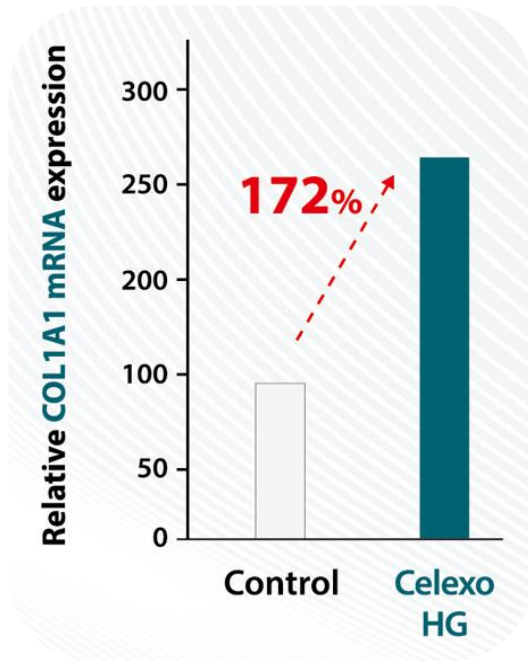

Collagen synthesis: ↑172%

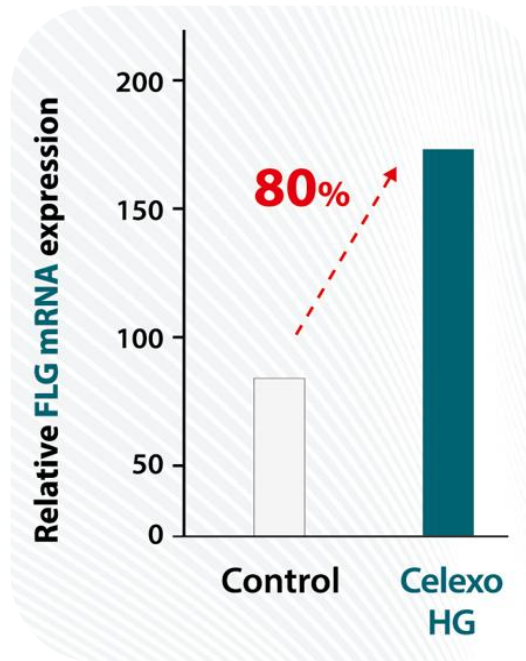

Skin barrier repair: ↑80%

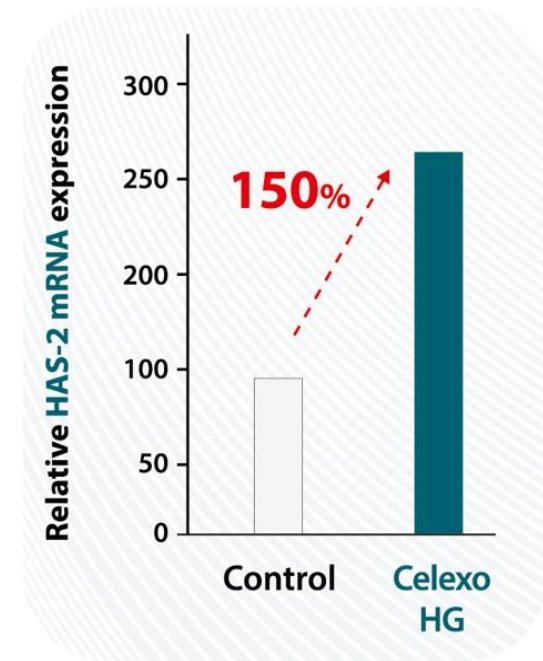

Hydration: ↑150%

Protocols

Step 1	Cleansing	Make sure the skin is completely free of impurities, makeup, and excess oils before starting the treatment.
Step 2	Anesthetic Cream	Apply an anesthetic cream for 20–30 minutes to ensure a comfortable treatment experience. Afterward, cleanse and disinfect the skin thoroughly.
Step 3	Select an appropriate device	Choose an EBD (Energy-Based Device) such as a laser or microneedle RF, based on the patient's skin type and condition.
Step 4	Celexo HG	Evenly apply an adequate amount of Celexo HG to the treatment area.
Step 5	Absorption Device	Use an absorption device to help the product penetrate deeply into the skin for maximum effectiveness.
Step 6	Soothing care	Gently apply a soothing face mask for 15–20 minutes to calm and relax the skin after treatment.
Step 7	Final skin care	Apply a regenerating cream and sunscreen to protect and nourish the skin.

Recommended Treatment Cycle

- Intensive care:
Every 2 weeks, total of 3-5 sessions
- Maintenance: Repeat every 6-12 months


Post-Treatment Precautions

- Avoid heavy exercise and saunas for 48 hours
- Use sun protection when going out
- Avoid touching treated areas (minimize irritation)

Protocols with Device

Item	Protocol 1	Protocol 2	Protocol 3
	MTS + Celexo Hydro gel	Potenza + Celexo Hydro gel	Virtue RF + Celexo Hydro gel
Protocol	• Apply Celexo 3cc after MTS • LDM for soothing purposes	• Use Potenza CP25 Tip • Manual Mode / Level 5 / Mono Polar / 2MHz / 700mJ • Apply Celexo 3cc after stamping • 400–500 shots • Velvet mask pack + LDM for soothing	• Use Virtue RF 36-Pin Tip • Manual Mode_Deep RF / Level 3 / Pulse 500 ms / 2MHz / Cooling: ON • Apply Celexo 3cc after stamping • 400–500 shots
Recommended Depth	• Forehead: 1.0 mm • Eyes: 0.75 mm • Mid-lower face: 1.75 mm 	• Forehead: 1.5 mm • Eyes: 1.25 mm • Mid-lower face: 2.0 mm 	• Forehead: 1.2 mm • Eyes: 0.7 mm • Mid-lower face: 1.5 mm 
	※ Treatment depth can be adjusted according to skin thickness		
Expected Results	• Anti-inflammatory effect for troubled skin • Accelerates skin barrier recovery • Provides soothing and calming effects	• Enhances exosome absorption • Improves acne & inflammation • Improves skin elasticity	• Stimulates collagen regeneration • Improves wrinkles & skin texture

Clinical Case | Sagging Skin

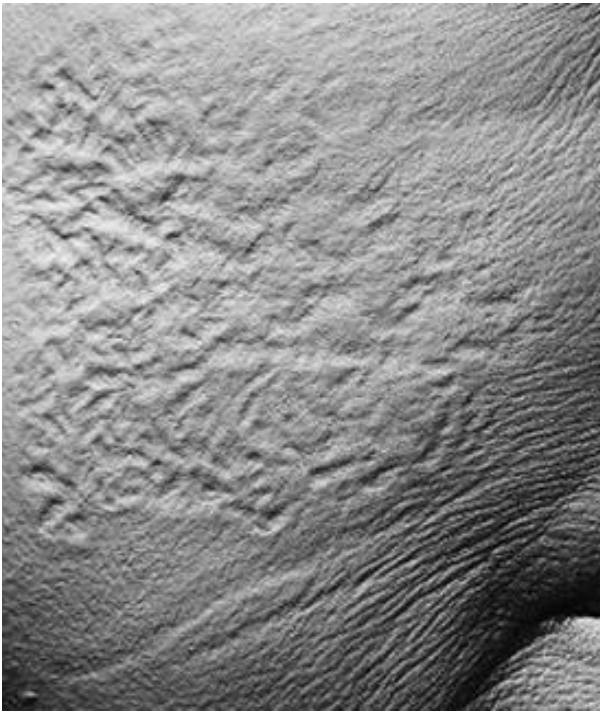

Before



After 3rd Session

Combined device – MTS

*Pillow friction test

Indication	Sagging Skin
Patient	Asian, Fmale, 40s
Protocol	Celexo Hydro gel 2ml with MTS
Results	Reduced roughness, faster recovery via collagen boost

Clinical Case | ACNE & Redness



Before

After 5th Session

Combined device – MTS

Indication	Acne, Redness
Patient	Asian, Male, 22
Protocol	Celexo Hydro gel 2ml with MTS
Results	Significant reduction in active acne lesions redness reduction

Clinical Case | **ACNE & PIH**



Before



After 5th Session

Combined device – **MTS**

Indication	Acne, PIH
Patient	Asian, Female, 21
Protocol	Celexo Hydro gel 2ml with MTS
Results	Significant reduction in active acne lesions Noticeable decrease in overall redness and PIH

Clinical Case | **Acne Scar, Redness**



Before

After 5th Session

Combined device – **Needle RF**

Indication	Acne Scar, Redness
Patient	Asian, Male, 45
Protocol	Celexo Hydro gel 2ml with MTS
Results	Improved skin texture, reduction in acne scar Noticeable decrease in redness

ABio materials

Advanced Bio materials

(주)에이바이오머티리얼즈

05854 서울특별시 송파구 송파대로 201, 제A동 1306호(문정동, 송파테라타워2)

Phone +82.2.6959.3734 | **E-mail** master@a-bio.co.kr

Copyright © 2025 Celexo. All Rights Reserved.

ABio materials USA

6262 Katella Ave. 2nd Floor Cypress, CA 90630

▼

www.a-bio.co.kr

By Exosome Professionals.

www.exosomeskinbooster.co.kr